




ONE-HANDED OPERATION INSTRUCTIONS

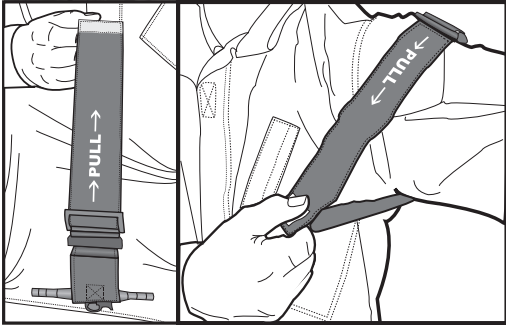


1 Hold the end of the strap by the white time label tab, allow the TMT to fall open with the word "PULL" facing away from you. Place around the injured extremity high and tight above the wound.

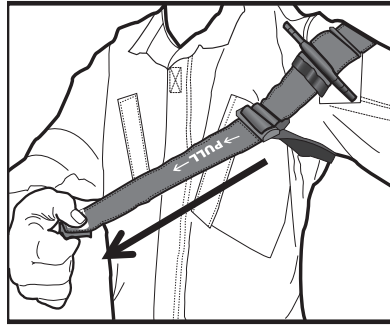


2 Pull time tab until Velcro® Brand material tail is separated from the TMT.

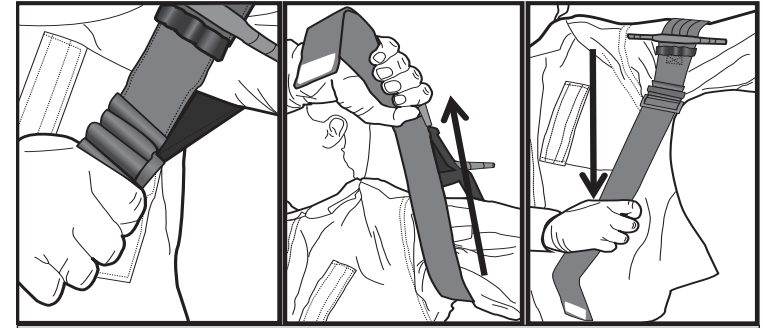


3 Grip tail close to the buckle, firmly pull up, out, and away from the body and then downward in a quick ratcheting motion. Repeat step until snug around extremity.

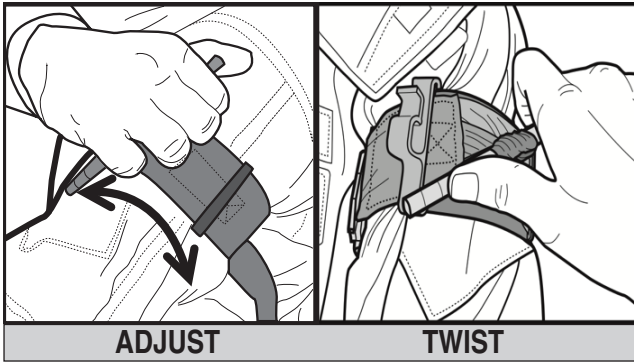

4 Adjust/rotate the TMT on extremity for ease of access to the torsion bar. Twist torsion bar until bleeding stops and the pulse is not present.

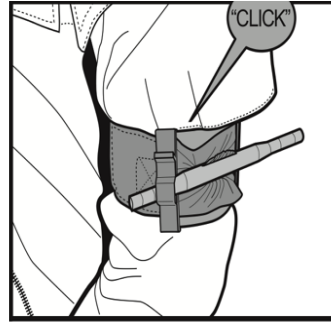

5 Lock the torsion bar into the retainer. Secure excess strap with Velcro® Brand material.

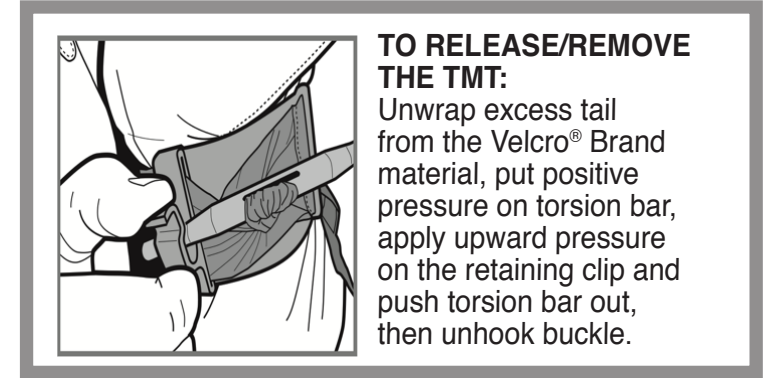

REMOVAL



TACTICAL MECHANICAL TOURNIQUET

NSN: 6515-01-656-6191
Cage: 056F9 PN# 60064



Always store in one-handed configuration.



Savior Global
8 Ash Gardens
Hamworthy, Poole, Dorset
BH15 4DE United Kingdom

FDA Registered
Patent Pending

This is a SINGLE USE device.
ISO 9001: 2008 Certified



Limitation of liability: Alphapointe, its employees, agents, contractors, suppliers, and distributors shall assume no liability for injury or damages arising from the application and use of the TMT. The user assumes all risk liability. The TMT should only be used as directed in the military service component guidelines, EMS authority, or under direct supervision of a physician.

Standard Use Guidelines: For life threatening hemorrhage situations. For the military in the Care Under Fire situation, treatment normally consists of using a tourniquet immediately to stop major bleeding of the extremities. In field care situation, the trained and authorized rescuer will perform all treatment necessary using the TMT as a last resort and should only be applied when bleeding cannot be stopped and the situation is life threatening.


Alphapointe
3611 14th Avenue
Brooklyn, NY 11218

Made in U.S.A.


TWO-HANDED OPERATION INSTRUCTIONS

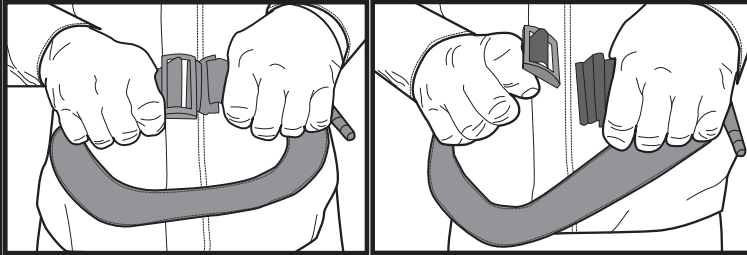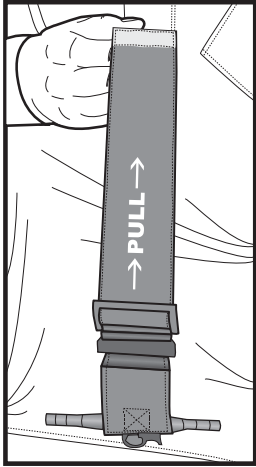


1 Hold the end of the strap by the white time label tab, allow the TMT to fall open. Grab TMT by the buckle and separate the two sections.

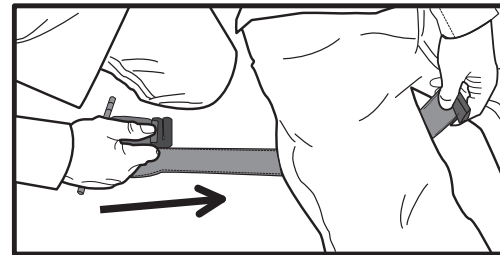


2 Wrap/Feed the TMT around the injured extremity high and tight above the wound, ensuring device is not twisted. Reconnect the buckle.

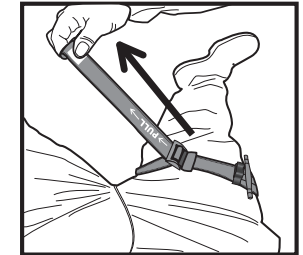

3 Pull time tab until VELCRO® Brand material is separated.

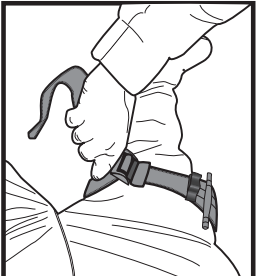

GRIP

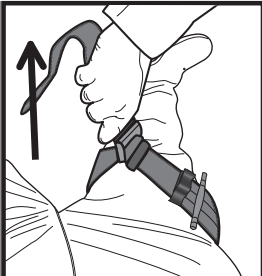


PULL UP




PULL DOWN

4 Grip tail close to the buckle, firmly pull up, out and away from the body and then downward in a quick ratcheting motion. Repeat step until snug around extremity.

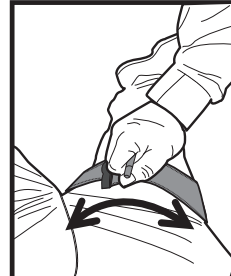


ADJUST

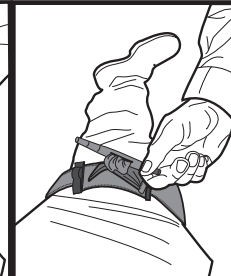

TWIST

5 Adjust/rotate the TMT on extremity for ease of access to the torsion bar. Twist torsion bar until bleeding stops and the pulse is not present.

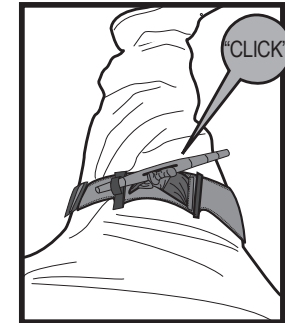

6 Lock the torsion bar into the retainer. Secure excess strap with VELCRO® Brand material.



**TACTICAL
MECHANICAL
TOURNIQUET**


NSN: 6515-01-656-6191
Cage: 056F9 PN# 60064



Alphapointe
3611 14th Avenue
Brooklyn, NY 11218



Made in U.S.A.

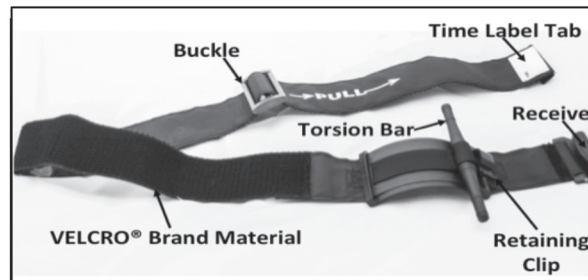

FDA Registered
Patent Pending

This is a SINGLE USE device
ISO 9001:2008 Certified

9 Line MEDEVAC Request

- Line 1: Location of the pick-up site
- Line 2: Radio frequency, call sign, and suffix
- Line 3: Number of patients by precedence
- Line 4: Special equipment required
- Line 5: Number of patients
- Line 6: Security at pick-up site
- Line 7: Method of marking pick-up site
- Line 8: Patient nationality and status
- Line 9: NBC Contamination (Peacetime: Terrain)